



Fisher Hamilton®



Tradition™...The beauty of being tough.

A lab is a very hostile environment for furniture. It is filled with spills, corrosive agents, heavy use and abuse of every description. To survive and last, lab furniture has to be built tough enough to take the worst that lab use can dish out. Steel Casework has become the product of choice for industrial labs, because it stands up to the challenge of years of intensive lab usage and heavy workloads while remaining beautiful and stylish. Fisher Hamilton Tradition Inset Steel proves conclusively that toughness and beauty are not mutually exclusive terms.





Create the perfect combination with Contrast™

Now it's possible to combine the rugged durability of Fisher Hamilton Steel cabinet construction with Wood or Plastic Laminate door and drawer fronts – all in contemporary Inset styling. Take advantage of the rich warm tones of Wood or the many color options of Plastic Laminate with all of the benefits of durable Steel furniture construction.



UL1805 Classified Fume Hoods

All Fisher Hamilton fume hood styles and sizes are performance tested to the UL 1805 standard.



Wall Cabinet Options

Wall cabinet options include open fronts, and sliding or hinged doors in solid panel, framed and unframed glass styles.



Work Surfaces & Fixtures

Choose from five different surfaces, and a wide selection of heavy-duty service fixtures to meet the specialized needs of your lab.



Durable Hinges

Durable, corrosion resistant stainless steel hinges.



Remote Leveler Access

Unique design permits leveler access through toe space area – no leveler access holes on cabinet floor pan.

Easy-to-Clean Drawers

Coved horizontal corners make drawers easier to keep clean.



Dependable Drawer Stop Design

Drawer stops contact welded drawer body rather than drawer head to provide longer, more trouble-free service.

Creative lab design is easy with modular components

Fisher Hamilton Inset Steel laboratory casework is available in a vast array of cabinet configurations and door and drawer front materials, all in Inset styling. Choose from Tradition or Contrast door and drawer front material options in an assortment of base, wall and upper cabinets, tall cases, specialty cabinets and utility tables. Then, complete the lab with work surfaces, sinks, fume hoods, fixtures and accessories.

Choice of Pull Styles

Rectangular brushed aluminum pulls are standard. Optional pulls in a variety of styles are also available.



Adjustable Shelf Height

Shelf height is adjustable in half-inch increments over the full height of the cabinet for the most flexibility in the industry.



Sleek Styling

Fisher Hamilton's smooth front narrowline styling features the narrowest stiles and rails in the industry to create a sleek, modern appearance.

Quality sets us apart



Wrap around cabinet body construction

The sides, back and outside cornerposts of each cabinet are made from one sheet of 18-gauge cold rolled steel to provide maximum strength and eliminate the need for an applied finished back.



Double corner posts at all four corners

Double corner posts are used to further enhance the strength and durability of each Fisher Hamilton cabinet. They provide the load bearing capacity to handle the heaviest laboratory equipment.



Ball bearing roller drawer suspension

Drawers operate almost effortlessly, and return to the same location every time for uniformity of reveal spaces. A lock open feature permits working from an open drawer with both hands free. Rise up drawer suspension prevents interference of a fully loaded drawer with an open door below.




Electrostatically applied urethane powder paint

Fisher Hamilton's SEFA 8 tested urethane powder coating provides reliable chemical and abrasion resistance. Select from 18 standard colors.

Door and Drawer Front Material Options

Steel Options

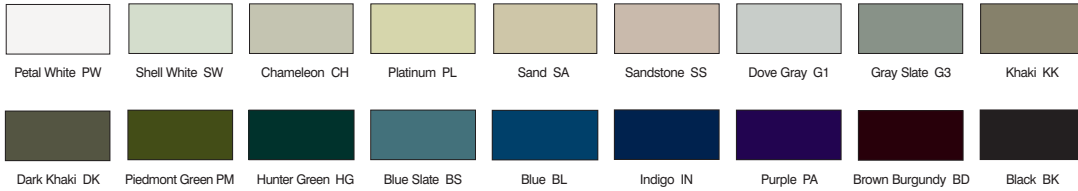

Paint Colors

Choose from 18 standard colors.

Wood Options



Stains/Finishes

Choose from seven finish colors.

*Not available in Hard Maple.

Veneers

Plastic Laminate Options



Colors

The 27 colors shown above are representative of the many Wilson Art plastic laminate colors in the Fisher Hamilton standard offering. Please consult with your Fisher Hamilton representative for availability.

NOTE: All colors are representative only. They are close representations of the actual colors within the limitations of printing. Actual samples are recommended prior to making a final selection.



PVC Edge Banding

Black, Gray, Chameleon, White

Common Options

Pulls



Rectangular Brushed Aluminum Finger Pull (Standard)

Wire Finger Pulls; Chameleon, Black, Gray, Brite Chrome and Goldtone Powdercoat; Satin Chrome, Polished Chrome Plated, Stainless Steel and Aluminum

Semi-Recessed Polypropylene Pulls; Black, Chameleon, Gray and White Aluminum Flush Pull

Full-Width Powder Coated Steel Flush Pull

Full-Width Wood Flush (Available on Wood Door and Drawer Fronts Only)

Full-Width Aluminum Flush Pull

Door Hinges



Black with Zinc-Plated Screws

Stainless Steel with Stainless Steel Screws

Drawer Suspension

SEFA 8-Self-Closing, self-centering, hold-open feature 100 pound drawer runs

Self-Closing, self-centering, hold-open feature 100 pound drawer runs – (Standard)

SEFA 8-Accuride 150 pound full-extension drawer slides (Option) (Standard on File Drawers)

Locks

National Master 5-Disc (Standard)

National Grand Master 5-Pin

Label Holders

None (Standard)

One Each Door, One Each Drawer Head; Black, Chrome or Goldtone

One per Pair of Doors, One Each Door Head; Black, Chrome or Goldtone

Removable Back Panel

Available on 100S and 200S Series Base Cabinets Only

At Cupboard Areas Only (Standard)

At Both Cupboard and Drawer Areas

Glass

Clear Float Glass (Standard)

Tempered Safety Glass

Laminated Safety Glass

Accessories

There is no need to go anywhere else when designing a lab furniture system. Fisher Hamilton provides all the tools needed to build the total lab environment. Here's just a sampling of the available accessories.



Work Surfaces Choose from five different work surface materials in a variety of colors, edge configurations and durability ratings.



Sinks Fisher Hamilton offers drop-in and self-rimming styles in a variety of configurations and materials. Epoxy resin is highly resistant to acids, solvents and salts. Stainless steel sinks feature polished, sound-deadened interiors.



Fixtures Choose from a variety of WaterSaver™ and BROEN® water faucets and service fixtures. WaterSaver fixtures are available in traditional chrome-plated or stylish ColorTech™ white powder coat finishes. BROEN fixtures feature a highly durable white powder enamel finish.



Fume Hoods The most technologically advanced and diverse line of fume hoods available today. From bench-type to floor mounted styles, Fisher Hamilton's high-performance fume hoods meet a full range of safety and ventilation needs.

Solutions for the Technical Environment

When it comes to choosing laboratory furniture, look to Fisher Hamilton. As the industry leader, Fisher Hamilton offers the most comprehensive selection of product lines, options and accessories in the world today.

The markets served are as diverse as the products manufactured. We work very closely with the educational, industrial and health care markets and are actively engaged in the development of new technologies and products to better serve these ever-changing laboratory environments.

The many lab furniture lines that can be drawn upon to build the optimal system include:

Wood Casework

Steel Casework

Plastic Laminate Casework

Phenolic Resin Casework

Stainless Steel Casework

MAX/Lab® Adaptable Furniture System

MAX/Wall® Technical Wall System

MAX/Mobile® Adaptable Furniture System

Distinction® Laboratory Bench System

C-Frame™ Adaptable Furniture System

P-Frame™ Adaptable Furniture System

System XL™ Adaptable System

Adaptable Table Systems

Laboratory Fume Hoods

BioSafety Cabinets

Dealer Planning Expertise

Fisher Hamilton's dealers are a network of full support laboratory specialists who prefer to get involved in the earliest planning stages possible so they can offer helpful suggestions that can make your laboratory work easier and more efficient than ever before.

Whether it's product renovation or entirely new construction, Fisher Hamilton dealers have the skills, the expertise and the product selection to get the job done efficiently and cost-effectively. Give your nearest Fisher Hamilton dealer a call today for assistance with the planning process.



Fisher Hamilton L.L.C.
1316 18th Street, P.O. Box 137
Two Rivers, WI 54241
Phone 920-793-1121
Fax 920-793-3084
www.fisherhamilton.com

